

FLEX

A POWERFUL AND FLEXIBLE SPECTROMETER THAT FITS IN THE PALM OF YOUR HAND

- *Perfect for General Purpose Applications*
- *Improved Optical Design – The Best in Class!*
- *User-Interchangeable Slits and Filters*
- *UV-Vis-NIR (185-1100 nm) Detection Range*
- *Fiber Optic Spectrometer*
- *Software LightScan Included*



PERFECT AS A GENERAL-PURPOSE SOLUTION

With a compact and flexible design, FLEX is more than capable of delivering a great performance with improved sensitivity and with the ability to be integrated in most general-purpose spectroscopic setups.



EASY TO INTEGRATE IN YOUR SETUP

The SMA fiber optic connector allows for a simple attachment of FLEX to a wide range of spectroscopic setups without the need for complex optical alignments. This flexibility is complemented with a trigger integration tool, which includes triggers in/out and delays. FLEX and all spectroscopic components of your setup can be centralized for full control with our new SPEC Hub solution.







CONTROL SENSITIVITY AND RESOLUTION OF YOUR SPECTROMETER

With the FLEX spectrometer, you can take advantage of interchangeable slits and filters. User-interchangeable slits offer an easy and simple way of setting the optimal resolution/sensitivity while interchangeable filters provide the user with the ability to narrow the operational wavelength range of the spectrometer. These features allow you to change from absorbance to a fluorescence setup in a few minutes. Add a SLIT KIT to your order and get all flexibility you need.

BEST OPTICS FOR A LOW STRAY LIGHT

FLEX is the perfect solution for those who have the need for absorbance/transmittance measurements. The simultaneous optimization of the optical and mechanical design provides FLEX with a very low stray light percentage (less than 0.1%, over the entire operational spectral range).

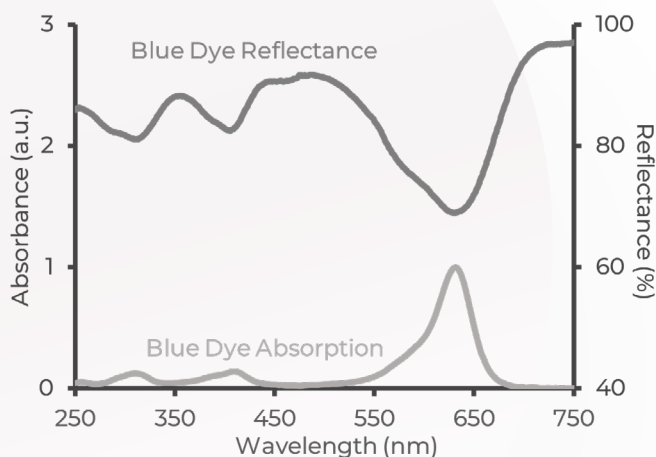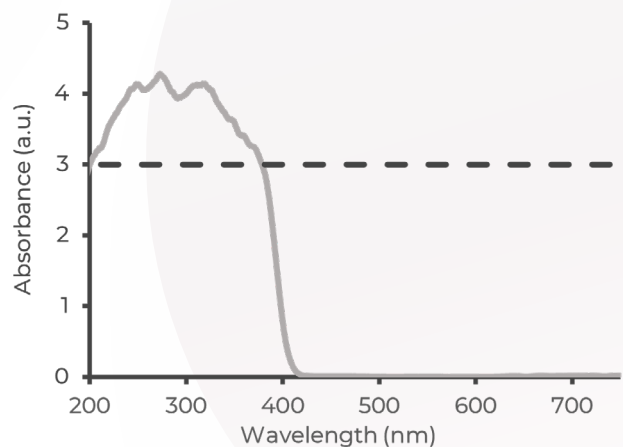

ABSORBANCE AND REFLECTANCE IN ONE SPECTROMETER

The versatility of FLEX can be easily achieved when combined with high-power deuterium and tungsten/halogen (DWHP) light source. In addition to an absorbance configuration (it requires the use of standard cuvette holder and optical fibers), this combination can also be used to perform reflectance measurements (it requires the use of reflectance probe).



+351 223 709 107



www.sarspec.com




info@sarspec.com



MAXIMUM FLEXIBILITY WITH THE BEST PERFORMANCE



Absorbance

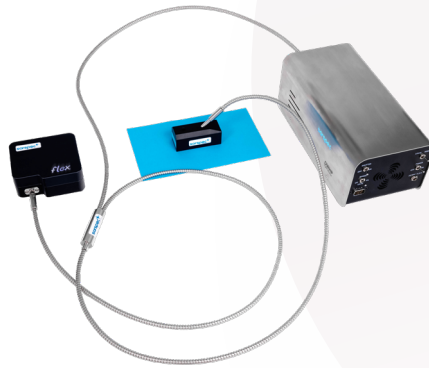

Reflectance



Fluorescence

FLEX is available in both standard (STD) and improved-resolution (RES+) versions for ready-to-use (table below) and user-configured configurations. For both solutions, users have the possibility of choosing an improved sensitivity configuration, which translates into a more efficient light collection optics. Also, customized configurations provides you with the flexibility to configure FLEX for more specific applications.

Technical Specifications

Parameter	UV-Vis-NIR	UV-Vis	Vis-NIR
Operational Range	185-1100 nm	185-900 nm	360-1100 nm
Grating	G-UV-Vis-NIR	G-UV-Vis	G-Vis-NIR
Op. Resolution STD	1.7 nm (25 um slit)	1.4 nm (25 um slit)	1.4 nm (25 um slit)
Op. Resolution RES	0.75 nm (10 um slit)	0.6 nm (25 um slit)	0.6 nm (25 um slit)
Fiber Connector	SMA-905		
Interface	Mini-UBS (works also as power supply))		
Integration Time	2 ms STD / 3 ms RES+ (minimum)		
Dimensions	192 mm (L) x 137 mm (D) x 75 mm (H) (Weight: 450 g)		

Ordering Information

Part Number STD	FLEX-STD-UV-Vis-NIR	FLEX-STD-UV-Vis	FLEX-STD-Vis-NIR
Part Number RES	FLEX-RES-UV-Vis-NIR	FLEX-RES-UV-Vis	FLEX-RES-Vis-NIR



+351 223 709 107



www.sarspec.com



info@sarspec.com

